MOVIX DREAMY





OVERVIEW



Movix DReamy: access the next generation of STEPHANIX digital mobile radiography.



Movix DReamy is designed to fulfil today's challenges and users' needs:

- ➤ Compact
- > Telescopic Column
- > Powerful (up to 50 kW)
- > Fully motorized
- Easy moving

x-ray.








OVERVIEW

Effortless handly thanks to the capacitive technology, user can easily adjust, move and position the mobile

Movix DReamy, a mobile designed for





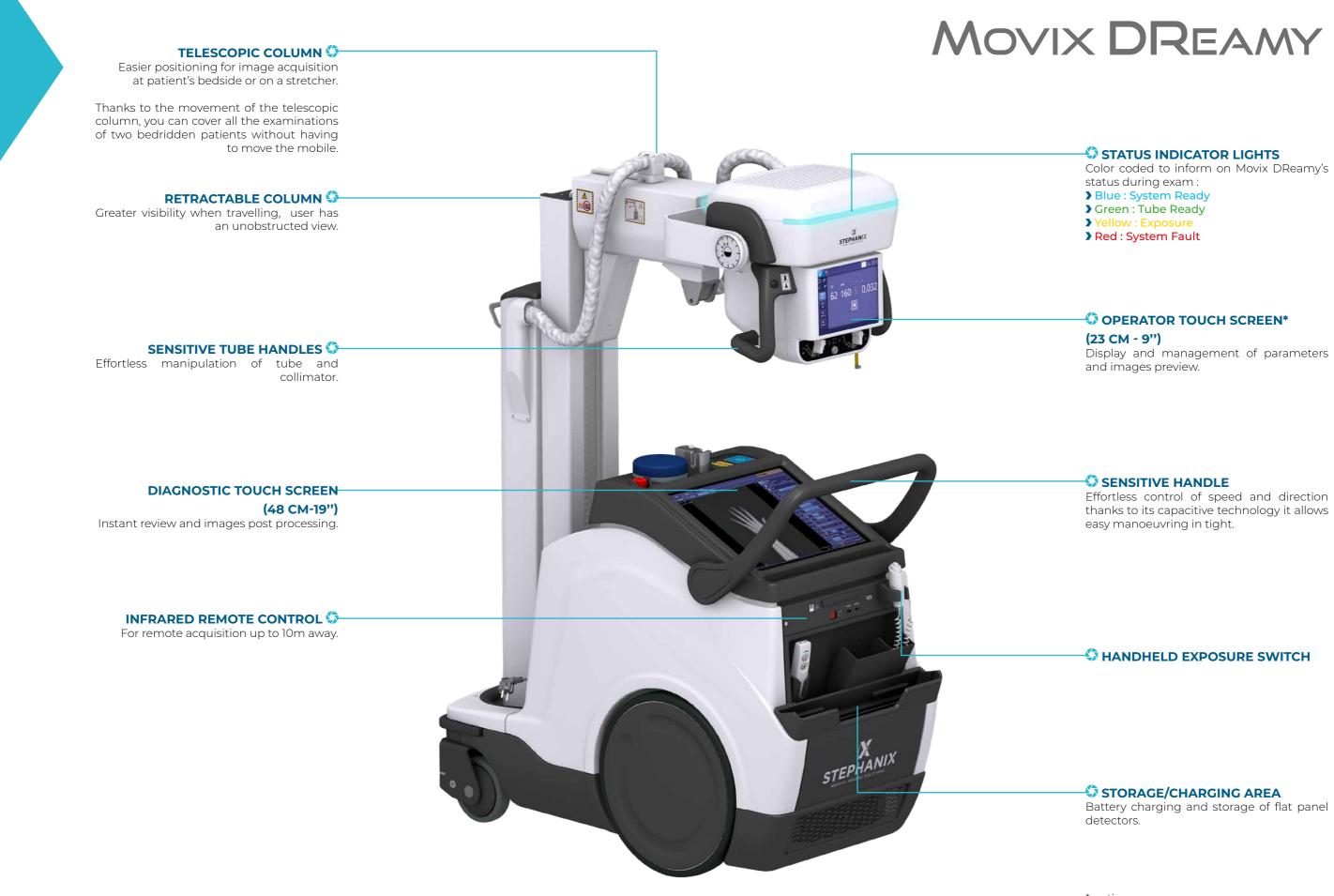

thanks to its capacitive technology it allows

* option

FREEDOM OF MOVEMENTS

Retractable column with a height of 129 cm, offers perfect visibility while moving.

Its compactness, with a width of 54 cm, guarantees maneuverability and easy access to the various rooms and departments.

Healthcare professionals can move around comfortably and safely in every corridor and care unit.









MOBILE MOVEMENTS

Moving speed is adjusted according to the pressure exerted on the handle and can reach up 5 km/h.

Move commands on tube head make it easy to adjust positioning while staying close to the patient. No need to use the drive handle for minor adjustments in close confines.

TUBE MOVEMENTS

Positioning in One movement: simultaneous combinations of rotations allow orientation of the tube and collimator in all axes.

Column, tube and collimator rotations provide best positioning for various types of exams.



Tube Horizontal Rotation +90/-30°



Tube Vertical Rotation : +/- 180° Mobile's movements on the tube head: Easier and faster positioning at patient's side without the need to return to the drive handle.



* Depending on configuration.

MOVIX DREAMY







Collimator Rotation : +/-90°

Movix DReamy rear diameter load wheels, as well as front steering wheels featuring shock absorbing springs help tp cross obstacles such as elevator threshold.

Easy crossing of obstacles, such as an elevator threshold

SAFETY FIRST

Movix DReamy has specific features to answer real life situation:

> Cramped corridors?

> Small and challenging doorways?

Tight rooms?

Specific Focal length?

> Radiation protection?







LED light field and focus laser

THE RIGHT S.O.D

By using LED light field and focused 1 meter double laser*, you can easily respect the distance between the x-ray source and the patient (S.O.D).

Guarantee patients and care givers maximum efficiency for a minimum dose, as well as compliance with radiation protection.

* Option, depending on configuration.



SECURE ACCESS BY BADGE (OPTION)

If user traceability and while allowing the charging of



TFT SCRREN* (OPTION) The on-tube display gives you acquisition quality.



INFRARED REMOTE CONTROL (OPTION)

turn on the light field and trigger reducing the scatter received by the

FEATURES









BUMPER

The safety bumper located on the mobile as soon as it collides with an obstacle.



STATUS INDICATORS LIGHTS Color Coded Leds indicate the status of your Movix DReamy.



DOSE MANAGEMENT Transfer of image dose information object to RIS/PACS compatible.

EASY OF USE

4

System start-up: Switch on the user interface. > Access worklist or manual entry.

BE MOBILE, STAY CONNECTED

Movix DReamy has an integrated console with a touch interface: you can select a patient from the worklist or perform a manual entry and choose one or more protocols.

The image acquisition and processing parameters are automatic and adapted to each protocol.

After exposure, visualization of the image is effective in less than 2s.

Various post-processing tools allow you to modify the contrast, add annotations, create measurements or even create a printing layout.

Scatter correction software will free yourself from the use of a physical grid. This will prevent misalignments and artefacts with the benefit of dose reduction.



Exam management: > Dose Display and Export > Images connectivity either via Wifi or

З

and user priority: > Dose and diaphragm management. > Manual filtration and scatter correction. > Optimized protocol.

DX MODALITY

SCATTER CORRECTION

Without correction





With correction










Movix DREAMY

2

Program

- Anatomical:
- > Presets Protocol Selection.
- > Automatic management of

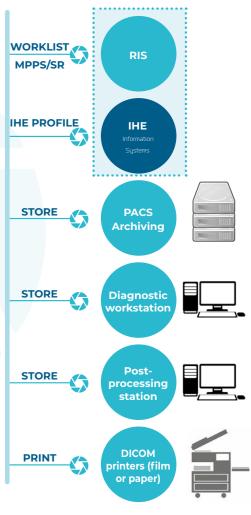

LATEST GENERATION OF DIGITAL RADIOGRAPHY SOLUTIONS: INTEGRATED, DYNAMIC, PORTABLE, WIRED AND WIRELESS DETECTORS.



X-RAY ROOMS

X-RAY MOBILE UNITS



SURGICAL C-ARMS





Head office: 10, rue Jean Moulin Z.I. du Bayon - 42150 La Ricamarie - FRANCE Tel. +33 4 77 47 81 60 - Fax: +33 4 77 37 55 19 export@stephanix.com

www.stephanix.com