

# RAD SERIES PRO DREAM



**STEPHANIX**  
MEDICAL IMAGING SOLUTIONS



Rad Series Pro DReam is a versatile RAD room to perform X-rays from head to toes. The exposures can be made directly on the table, but also on a stretcher or wheelchair.



The version with variable height adjusts the access of the tabletop to the stretcher for a better patient transfer.



**COMFORT OF USE**

The four-ways floating tabletop is made of low absorption materials < 1.0 mm Al eq.

Its width is suitable for patient positioning up to 300 kg.



The Wall Bucky moves 150 cm vertically to explore the entire skeleton.

In auto-tracking configuration, the tube movements are synchronised with the vertical displacement of the Wall Bucky, and with the table whenever applicable.


A tilting Bucky is available for extremities exams.

## POSITIONING

Both vertical & horizontal exposures are possible thanks to the movements of the column tubestand. Travel on rail also offers a wide X-ray coverage.

The +/- 180° column rotation allows direct projections on stretcher: the patient transfer to the table is no longer necessary.

The x-ray assembly rotations, controlled by electromagnetic brakes, enable the positioning to realise many types of exams.

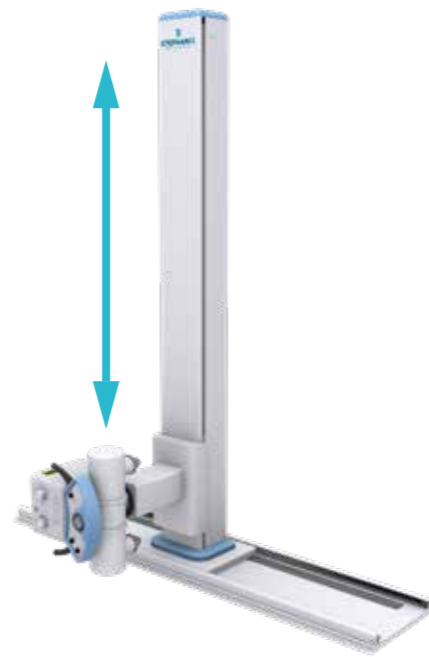

Rad Series Pro DREAM is equipped with a touch screen in auto-positioning configuration on which the acquisition parameters can be modified.

The functions managed by auto-positioning are:

- ▶ Angulation
- ▶ Movements
- ▶ Collimation and automatic filters
- ▶ Automatic centering
- ▶ Automatic focal distance



### FOCAL DISTANCE

Adjustable SID with displacement on rail to perform chest exams.



### VERTICAL MOVEMENT

Vertical travel of the Wall Bucky and tube for exams like lateral weight bearing foot and sinus in standing position.



### ROTATION

Rotation of the column tubestand on several axes for lateral projections.



### TELESCOPIC ARM

Telescopic arm for a useful transversal movement during examinations on polytrauma or on stretcher.



### TOUCHSCREEN (OPTION)

Touch screen in auto-positioning configuration with display and modification of the acquisition parameters available.



### DIGITAL DR

Flat panel technology with real time control of image quality.



### FLAT PANEL DETECTORS

The flat panel detectors, integrated wired or wireless, are designed for exams such as chest, extremities, abdomen and spine.

WiFi detectors thickness is equivalent to CR cassettes and enables specific procedures such as femoro-patellar and weight bearing feet.

Their sharing, between several rooms and mobile systems, permits all routine procedures.

### ACQUISITION

For digital versions with flat panel technology, a console with touch screen interface (option) gives access to the patient data entries, the selection of anatomical protocols and acquisition parameters, the image preview for validation, and its post-processing before sending it to the network.

For information sharing, **Rad Series Pro DReam** systems are equipped with:

- ▶ USB ports
- ▶ Ethernet outputs
- ▶ DICOM Media Export compatibility.




### WORKSTATIONS

Depending on your Department configuration, STEPHANIX is able to offer a selection of workstations for display and secondary processing with **Rad Series Pro DReam**.

STEPHANIX solutions manage the patient data, from image acquisition to diagnosis, and storage on the network.



## DX MODALITY

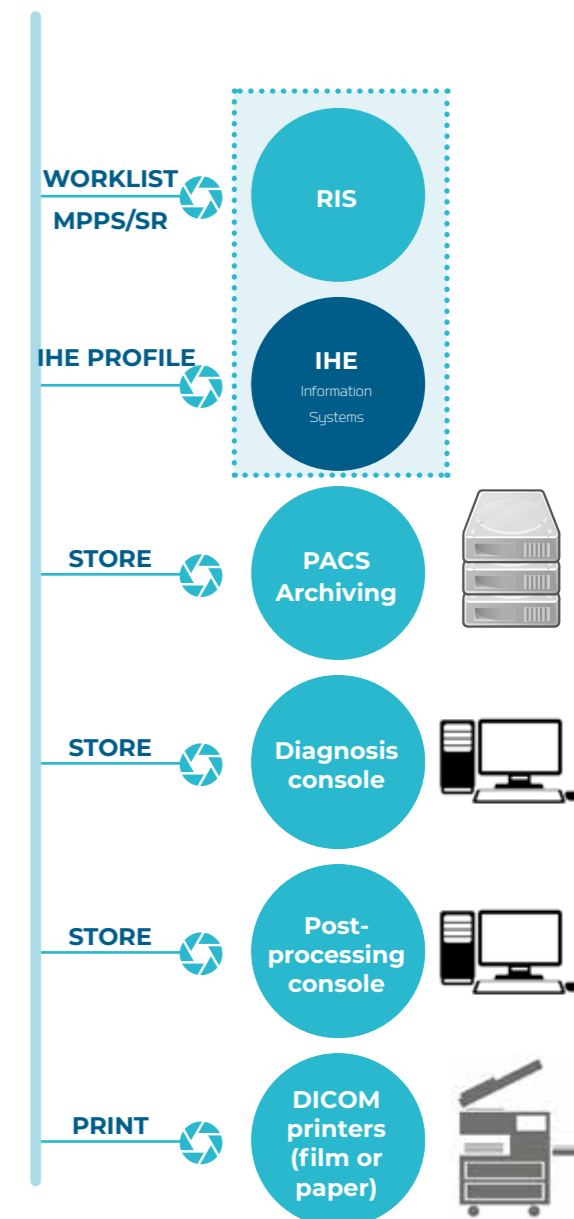

**DICOM**  
IT ARCHITECTURE GIVEN AS AN EXAMPLE

# LATEST GENERATION OF DIGITAL RADIOGRAPHY SOLUTIONS: INTEGRATED, DYNAMIC, PORTABLE, WIRED AND WIRELESS DETECTORS.

## REMOTE CONTROLLED TABLES



## X-RAY ROOMS



## X-RAY MOBILE UNITS



## SURGICAL C-ARMS



Real Series Pro Dream - 2022/05 - STEPHANIX se réserve le droit de modifier les caractéristiques de ses produits en fonction de l'évolution perpétuelle des technologies. Photographies non contractuelles. - Stephanix allows the right of modifying the products features regarding the evolution of the technologies without notice. Non-binding document.



**STEPHANIX**  
MEDICAL IMAGING SOLUTIONS

Head office: 10, rue Jean Moulin  
Z.I. du Bayon - 42150 La Ricamarie - FRANCE  
Tel. +33 4 77 47 81 60 - Fax: +33 4 77 37 55 19  
export@stephanix.com

[www.stephanix.com](http://www.stephanix.com)